

MANUAL IRRIGATION SET

The manual irrigation set is a medical device for transferring coolant (saline solution at 4-5 ° C) from the bottle to the patient's oral cavity. It is used with bag presses or with phisiodispenser.

DEVICE CLASSIFICATION

CLASS IIA	CE NUMBER: 0123	STERILE	Code: 270630
CND: A03010105 (national classification of medical devices)			
GMDN: Not applicable			

TECHNICAL FEATURES

- PVC tip to facilitate positioning on the cooling needle of the handpiece;
- Steel needle
- The product is free of latex and phthalate components (DEHP).
- The product is non pyrogenic.

Compatible with	Total Length	Pieces per pack
NSK SURGIC XT	226,4 m	10

APPLICABLE STANDARDS

ISO 8536-4, ISO 8536-8, ISO 10993-1, ISO 10993-7

PACKAGING

Irrigation sets in protective wrapping, individually packaged in bags/blisters for peel-open sterilization marked with the sterilization batch and packed in cardboard containers.

Keep away from sunlight. Keep dry.

STERILIZATION

Sterilization is carried out with ethylene oxide; the sterilization process is regularly validated and kept under constant control by routine analysis; and the residue of the sterilizing gas is kept within the limits established by the current legislation. Sterilization is performed by a specialized laboratory.

EXPIRY

The validity of the sterile device is **59 months** from the date of sterilization (MONTH-YEAR).
The expiry date indicated refers to the device in its unopened and correctly stored package.

LABELLING


Device with removable adhesive labels to facilitate traceability.
Each blister comes with removable adhesive labels to help operators with traceability of the devices used during operations. The removable parts of the label can be detached from the pack and applied to the medical record or to any other document that, according to the procedures used, is used for traceability.

In addition, each individual pack is labelled with all the information required to identify the product:

- | | |
|--|--|
| <ul style="list-style-type: none"> • product description (contents list) • trade name of the product • product code number • manufacturer's name • method of sterilization • production batch number | <ul style="list-style-type: none"> • expiry date • sterile symbol • single-use symbol • production area • CE marking and number of the notified body • Precautions |
|--|--|

There is an adhesive label on the transport packaging bearing all the information required to identify the product as indicated above and the no. of pieces per pack



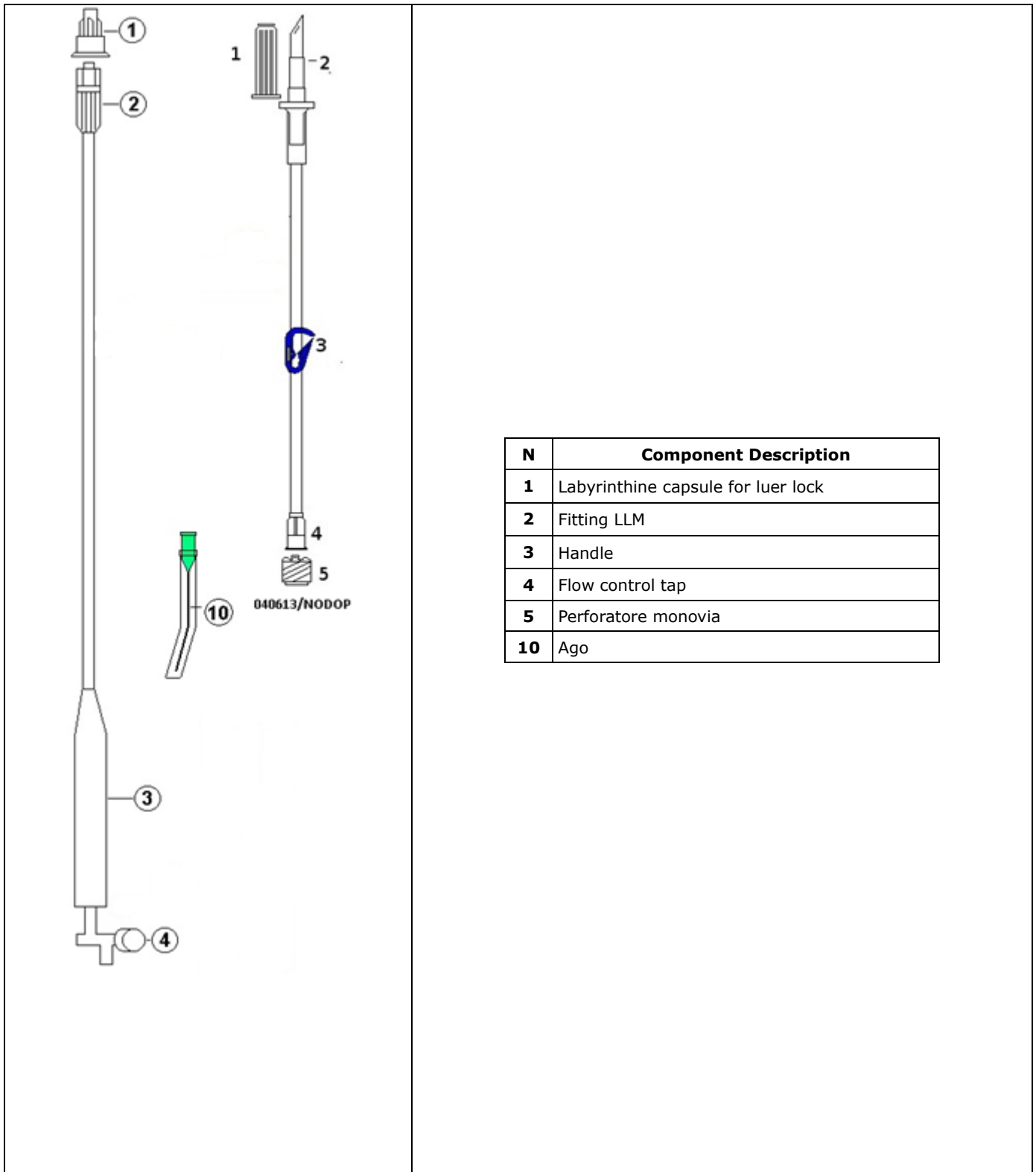

N	Component Description
1	Labyrinthine capsule for luer lock
2	Fitting LLM
3	Handle
4	Flow control tap
5	Perforatore monovia
10	Ago

